

ARJO

MAXI LITE

A PORTABLE LIFTING SOLUTION



ARJO'S NEW PORTABLE LIFT



In the emergency room, patient homes and nursing facilities, the *Maxi Lite* is the ideal lift for all daily lifting requirements. It performs general lifts, such as bed to chair, chair to toilet, chair to car, and from the floor, with ease.

The *Maxi Lite*, with a safe working load (SWL) of 350 lbs, is delivered to you assembled. The *Maxi Lite* has an integrated charger allowing for easy charging. For higher usage applications, there is an optional wall charger available, which assures the lift is always ready for use with a fully charged battery.

The slings are very easy to apply from the seated or supine position with the split leg design. They are safe and secure due to the no swing and sway design utilizing short clips instead of long straps. This allows the *Maxi Lite* to easily maneuver while the legs are in the closed position. The slings are also extremely comfortable with the padded leg sling design. This, together with the reclining action of the unique frame, creates a feeling of complete security. All of the slings are color-coded to make size selection easier.

The powered remote handset operation allows the operator to stay close to the resident at all times, which is a particularly useful feature when lifting from the floor. The *Maxi Lite* has completely dual controls for ease of use.

The *Maxi Lite* is collapsible, which means easier portability and less storage space. At its narrowest, the *Maxi Lite* passes through a standard doorway. At its widest, it fits around an armchair or wheelchair.

Features:

- Integrated charger
- Integrates with most auto types
- Optional Wall charger
- Adjustable legs
- Folding feature for storage
- Emergency shut-off power
- Front and/or side entry

Product information

Max. safe working load	350 lbs (160 kg)
Total lift height	Min.: 50.3" (1278 mm) / Max.: 71.1" (1805 mm)
Base total length	42.9" (1090 mm)
Base height	4.4" (112 mm)
Base floor clearance	1.9" (48 mm)
Battery	Two 12V, 5Ah (24V capacity)
Protection class	Handset IP44 Lifter IP21


Emergency stop and lowering

Approvals

The product is in compliance with current legislation and applicable product standards.

Resident Gallery

Suitable for Doris and Emma

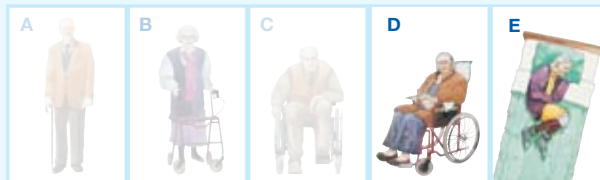

D Doris, who has no capacity to support herself

E Emma, who is almost completely bedridden and totally dependant

Please contact ARJO for further information on the Resident Gallery™.

Use only ARJO Designed Parts, i.e. parts specifically designed for the purpose, on equipment and products supplied by ARJO. Our Conditions of Sale specify that ARJO is not liable for injuries attributable to the use of inadequate parts. Continuous development is part of our policy, and we therefore reserve the right to change specifications without notice.

ARJO, *Maxi Lite* and Resident Gallery are trademarks belonging to the ARJO group.



...with people in mind.

ARJO, Inc.

50 N. Gary Avenue
Roselle, Illinois 60172
Phone: 800-323-1245
Fax: 888-594-2756
www.arjo.com

MEMBER OF THE GETINGE GROUP